



Legend Pro

Dare to imagine



Treat your skin

The complete platform for personalized aesthetic solutions that meet your patients' expectations and more

Imagine a treatment with immediate results, long-term efficacy, minimal pain and downtime, and no need for anesthesia...

Sounds too good to be true? **Legend Pro™ is all that, and more.**

The platform with the power of synergized technologies

Legend Pro™ offers endless options for personalized treatment plans to answer the specific needs and desires of each patient. With **Legend M2™** combined treatment, you can get optimal results and maximize patient satisfaction.

- Minimal pain and downtime
- Immediately visible results
- Proven combined therapy protocol
- Short treatment time



Make your patients' wishes come true while ensuring their safety and long-term satisfaction

Each patient is a world of their own; complete, unique and singular.

Legend Pro™ allows them to enjoy that individuality by enhancing their best traits with treatments bespoke to them.

The Legend Pro™ platform offers a complete tool kit for tailoring the right treatment to each of your patients, with the flexibility of treating a wide range of skin types, ethnicity and imperfections.

Treatment with **Legend M2™** reaches multiple layers of the skin, introducing a wide range of solutions for skin resurfacing and inducing renewal, as well as wrinkle reduction.

Testimonials

"The TriPollar® technology generates energy and ultimately remodels collagen for skin rejuvenation. The use of TriPollar® allows for heat to be disseminated through the skin with minimal pain and discomfort. Less pain means more comfort for our patients, and combined with the clinical efficacies which have been demonstrated, makes the Legend one of those go-to machines for skin rejuvenation. With the Legend system, I can totally achieve the results my patients and I desire. The facial dermal volumizing and wrinkle reduction results are exceptional. The skin becomes more supple and elastic and my patients look and feel years younger. Most importantly, the results are natural-looking which is what my patients want."

Michael H. Gold, M.D., Nashville, Tennessee, USA

"The TriPollar® RF and the VoluDerm™ allow me to customize the treatments for my patients and to get optimal results. My patients are happy going through almost painless treatment, and the treatment time is shortened without the need for numbing. These unique technologies allow me to treat multiple indications and get exceptional results."

Matteo Tretti Clementoni, MD – Laserplast Sagl, Milan, Italy



Courtesy of Steve D. Shapiro, M.D.



Before

After 8 treatments



Courtesy of Cruzira Salazar, M.D.



Before

After 2 treatments



Courtesy of Cruzira Salazar, M.D.



Before

After 3 treatments


By encouraging and enhancing the natural processes of the body, Legend Pro™ unveils the best in each patient.

From resurfacing of the upper layers of the skin to promoting the renewal of collagen and elastin fibers in the deeper dermis, Legend Pro™ offers a complete toolkit of non-invasive technologies for creating tailor-made treatments and personalized results for your patients.




TriPollar®

Triggering renewal of collagen and elastin fibers, by delivering focused RF energy in a safe and effective manner, inducing controlled heating of the deep dermis.

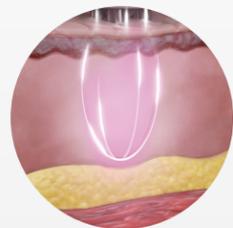

VoluDerm™

Activating the body's natural healing processes by applying RF assisted thermal ultra-thin electrodes ablation, creating micro channels to induce controlled micro wounds, thus encouraging the generation of new and healthy skin.

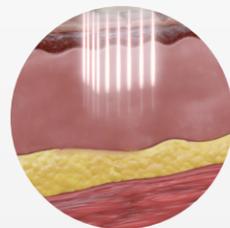

Smart Tip Design

Unique smart tip design enabling even distribution of RF energy that flows between the ultra-thin electrode and the tip surface, thus reducing the risk of arching and undesired side effects.




TriFractional™

Applying contact RF and pin point ablation with controlled volumetric heating to achieve skin resurfacing and renewal via wound healing process, resulting in skin texture improvement and refinement.



Built for perfection for both patients and physicians

Legend Pro™ includes unique design features that are created to ensure the safety and comfort the treatment to patients and physicians.

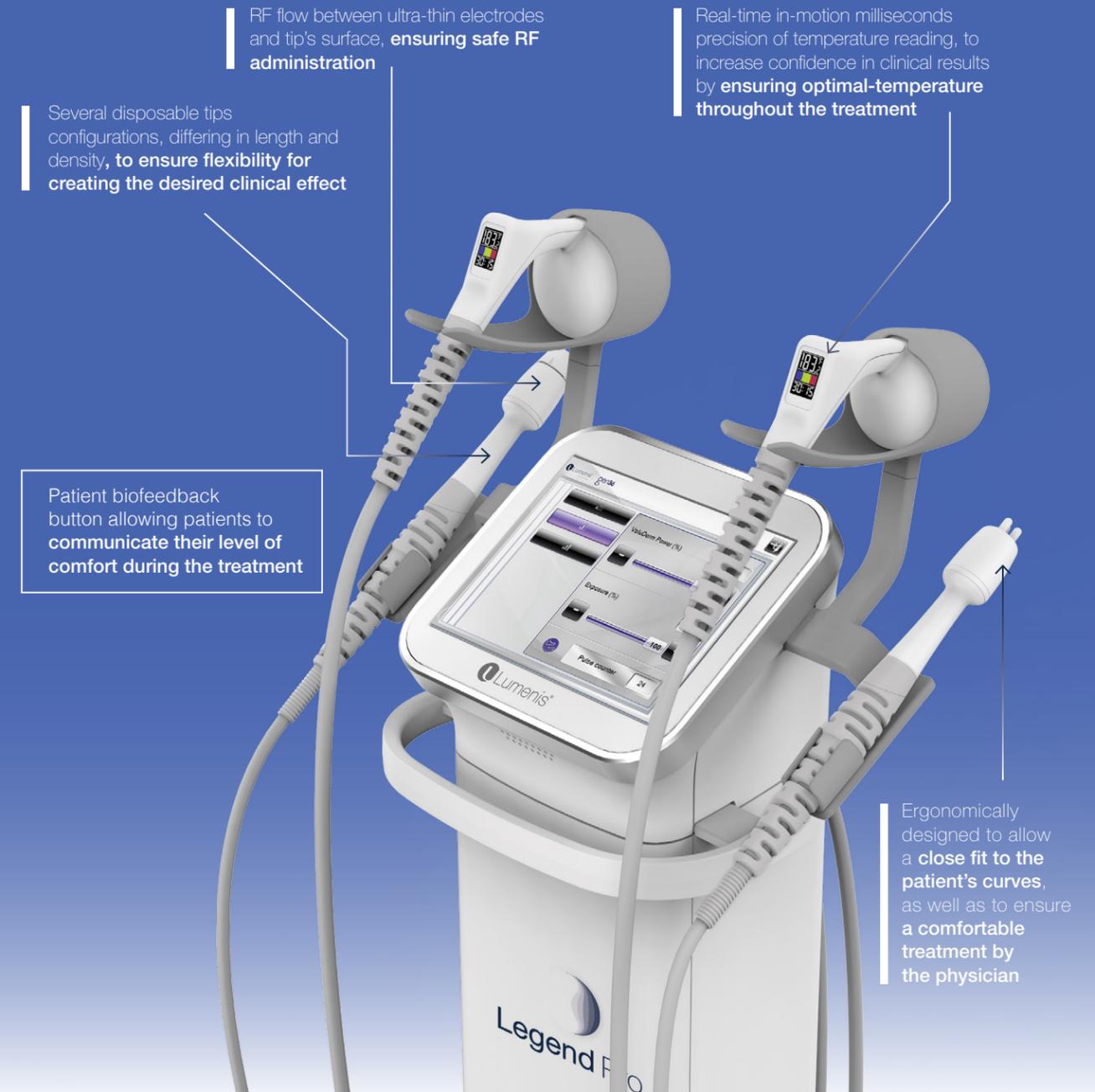

Legend M²™ clinically proven and tested combination therapy

The proprietary technologies in Legend Pro™ are scientifically proven through clinical research and market experience.

The TriPollar® RF and VoluDerm™ technologies target different skin layers individually. Preheating with the TriPollar® technology lowers the tissue impedance and increases vasodilatation, which may improve the efficacy of the VoluDerm™ ultra-thin electrodes.

During the Legend M²™ treatment the fibroblasts undergo two different triggers, minimally invasive, and non-invasive. The dual effect creates a synergetic result.

Legend M²™ unique synergized treatment combines the **macro** volumetric heating of tissue with TriPollar® RF and the precise **micro** heating of focal points with VoluDerm™ RF assisted ultra-thin electrodes to turn your patients' wishes into reality.

Seeing is Believing



Courtesy of Mark Rubin, M.D.



Before After 5 treatments

Courtesy of Mark Rubin, M.D.



Before After 3 treatments

Courtesy of Steve D. Shapiro, M.D.



Before After 3 treatments

Legend M²™ treatment validated through clinical experience, as published in the article by Steven D. Shapiro "Voluderm micro-needle Technology for Treating Skin Laxity and Wrinkles-Initial Clinical Experience".

Jacobs Journal of Experimental Dermatology, January 2015

Your journey with Lumenis has just begun.

We want your business to grow with Lumenis. That's why we have put together a support package that begins with the first conversation and will stay with you as your business begins to see the benefits. This ensures you have all the training you need and receive the most up to date clinical materials. Our direct to consumer marketing activities support the brand, keeping Lumenis Legend Pro™ the treatment of choice when it comes to wrinkle and resurfacing treatments.



Lumenis Practice Locator

Thousands of patients search for treatment providers every month. Help them find your clinic by registering your practice on the Lumenis Practice Locator. We run campaigns through influencers, social media and viral videos to promote this directory.



Lumenis online

Consumers have many questions about aesthetic. Our patient-facing websites will give them the re-assurance and confidence they need to move forward with the treatments.



Team Training

We install the system and train your staff, so that you are able to assess skin types and conditions, creating the perfect treatment plan for each individual patient. Our expert clinical training means they will handle consultations and treatments with confidence.



Your Most Valuable Resource

Your dedicated representative is an email, text, or phone call away. We want you to feel part of the Lumenis family, so please feel free to contact us any time. We'll always do our best to get back to you within 24 hours.



Brand Your Clinic

As a Lumenis client you will have access to the Lumenis Physician Portal. Here, you can find materials such as brochures and posters. You will also have access to a wealth of clinical and product information.



Extensive Service Network

We are committed to providing high quality, comprehensive service for medical equipment repair and maintenance. Lumenis service provides you with peace of mind, so you can focus on patient care, rather than on the equipment.

Specifications

Legend Pro™

Parameter	Value/Data
Input voltage	100-240V~, 50/60 Hz, max 2.2A
Mode of operation	VoluDerm™, TriFractional™, TriPollar®
Preset Parameters	Yes
Maximum output power RF	50 Watts @ 200 Ohm
Maximum output voltage	100±10 Vrms
Output frequency RF	1MHz
Platform Specifications	
Weight	~30 Kgs / 66 lbs
Dimensions	L18xW18 x H43 Inch

TriPollar® Applicator

TriPollar®	Number of electrodes	Temperature sensor
Applicator 1	6	Internal
Applicator 2	3	Internal
Applicator 3	3	External

VoluDerm™ Applicator

VoluDerm™ tips	Ultra-thin ectrode array	Length	Pulse count
Gen36	6x6	0.6 mm	Up to 800
Gen12	6x2	0.6 mm	Up to 800
Gen100	10x10	0.6 mm	Up to 800
Gen36L	6x6	1.0 mm	Up to 800
TriFractional™ tips			
	Pins array	Length	Pulse count
H7X7	7x7	0.2 mm	Up to 800

Warnings and risks:

Legend Pro™ is contraindicated for patients with pacemaker, defibrillator, or any implanted electronic device, metal implants in the treatment area. Improper use of the System could result in possible side effects.

Although these effects are rare and expected to be temporary, side effects may appear either at the time of treatment or shortly after. Side effects may include any of the following: prolonged or significant pain damage to natural skin texture (blister, burn), excessive skin redness (erythema), excessive swelling (edema), fragile skin bruising, excessive itching, change of pigmentation (hyper-pigmentation or hypo-pigmentation), scarring transient skin break-out such as acne and pimples. For complete list of contraindications, please refer to the User Manual.

About Lumenis®

Lumenis is a global leader in the field of minimally invasive clinical solutions for the surgical, ophthalmology and aesthetic markets, and is a world-renowned expert in developing and selling energy-based technologies, including laser, intense pulsed light (IPL) and radio frequency (RF).

For over 50 years, Lumenis' ground-breaking products have redefined medical treatments and have set numerous technological and clinical gold standards.

© 2018 All Rights Reserved. The Lumenis Group of Companies. Lumenis, its Logo, Legend Pro, VoluDerm, TriPollar and Trifractional are trademarks or registered trademarks of The Lumenis Group of Companies. Specifications are subject to change without notice. VoluDerm, TriPollar and Trifractional are proprietary technologies of The Lumenis Group of Companies.

Lumenis LTD

Yokneam Industrial Park
Hakidma Street 6,
P.O.B. # 240
Yokneam 2069204, ISRAEL
T: +972 4 959 9000

AMERICAS

San Jose, CA, USA
T: +1 408 764 3000, +1 877 586 3647
F: +1 408 764 3999



Pollogen Ltd.,

6 Kaufman St., P.O.Box #50320
Tel Aviv 6801298, Israel